




ImmunoAffinity Column for Zearalenone for sample preparation

Order no.: IC-CF-01-25
Lot: xxx xxx xxx xxx



aokinImmunoClean CF ZON column

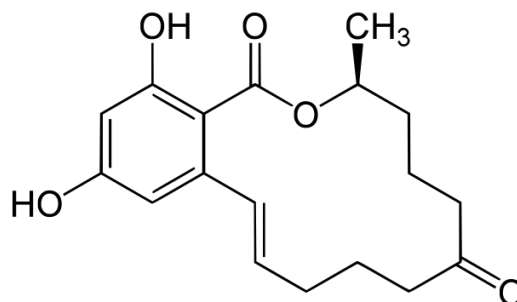

Zearalenone (ZON)

Specification:

Recommended loading: < 500 ng ZON
Loading rate: 1-2 ml/min
Elution rate: 1 ml/min

	Recovery of 0.2 µg
Specification:	<input checked="" type="checkbox"/> Recovery > 80%

Methodology:

Sample loading	ZON standard	0.2 µg Zearalenone
	Fluid	10% Methanol / 90% PBS, 0.05% Sodium azide, pH=7.4
	Volume	10 ml
	Rate	1 ml / min
Column washing:	Fluid	10% Methanol / 90% PBS, 0.05% Sodium azide, pH=7.4
	Volume	5 ml
	Rate	1 ml /min
Elution:	Eluant	100 % Methanol
	Volume	1.2 ml
	Incubation	3 min
	Rate	1 ml / min
Analysis:	by FPIA	



Dipl.-Ing. Biotechnology
Nicole Abdallah
Head of Quality Control