


## ImmunoAffinity Column for Ochratoxin A for sample preparation

Order no.: IC-CF-04-25  
Lot: xxx xxx xxx xxx



aokinImmunoClean CF OTA column

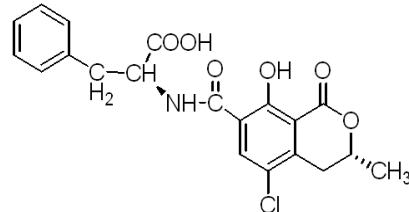

Ochratoxin A (OTA)

### Specification:

Recommended loading: < 200 ng OTA  
Loading rate: 0,5 - 1 ml/min  
Elution rate: < 1 ml/min

	Recovery of 65 ng	Recovery of 200 ng
Specification:	<input checked="" type="checkbox"/> Recovery > 80%	<input checked="" type="checkbox"/> Recovery > 80%

### Methodology Testing:

**Sample loading:** OTA 0 ng, 65 ng, 200 ng Ochratoxin A (OTA)  
Fluid 10% Methanol / 90% 10mM PBS pH 7,4  
Volume 10 ml  
Rate 1 ml / min

**Column washing:** Fluid 10% Methanol / 90% 10mM PBS pH 7,4  
Volume 10 ml  
Rate > 1 ml /min

**Elution:** Eluant 2% Glacial Acetic Acid / 98% Methanol  
Volume 1.5 ml  
Incubation 2 min  
Rate 1 ml / min

**Analysis:** by FPIA



Dipl.-Ing. Biotechnology  
Nicole Abdallah  
Head of Quality Control