




ImmunoAffinity Column for T2/HT2 toxin for sample preparation


Order no.: IC-C+-78-25
Lot: XXX XXX XXX XXX



aokinImmunoClean C+ T2/HT2 column

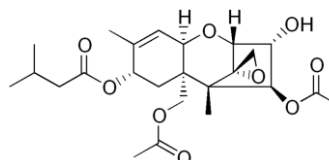

T2-Toxin (T2)

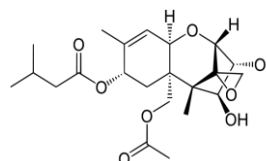

HT2-Toxin (HT2)

Specification:

Recommended loading: < 500 ng T2/HT2
Loading rate: 1-2 ml/min
Elution rate: 1 ml/min

	Recovery of 0,25 µg	Recovery of 1 µg
Specification:	<input checked="" type="checkbox"/> Recovery > 80%	<input checked="" type="checkbox"/> Recovery > 80%

Methodology Testing:

Sample loading:	T2/HT2 standard	0 ng, 0,25 µg, 1 µg T2/HT2 Mix 1:1 (T2, HT2)
	Fluid	10% Methanol / 10 mM PBS, pH 7.2 / 0,05% Sodiumazide
	Volume	25 ml
	Rate	1 ml / min
Column washing:	Fluid	10 mM PBS, pH 7.2 / 0,05% Sodiumazide
	Volume	5 ml
	Rate	> 1 ml /min
Elution:	Eluant	100% Methanol
	Volume	1.5 ml
	Incubation	2 min
	Rate	1 ml / min
Analysis:	by HPLC	



Dipl.-Ing. Biotechnology
Nicole Abdallah
Head of Quality Control