


ImmunoAffinity Column for Ochratoxin A for sample preparation

Order no.: IC-M-04-50
Lot: xxx xxx xxx xxx



aokinImmunoClean M **OTA** column

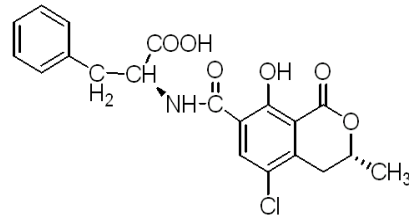

Ochratoxin A (OTA)

Specification:

Recommended loading: < 50 ng OTA
Loading rate: 1-2 ml/min
Elution rate: < 1 ml/min

	Recovery of 10 ng	Recovery of 200 ng
Specification:	<input checked="" type="checkbox"/> Recovery > 80%	<input checked="" type="checkbox"/> Recovery > 80%

Methodology Testing:

Sample loading: OTA 0 ng, 10 ng, 200 ng Ochratoxin A-Standard
Fluid 10 mM PBS, 10 % Methanol
Volume 10 ml
Rate 1 ml / min

Column washing: Fluid 10 mM PBS, 10 % Methanol
Volume 5 ml
Rate > 1 ml /min

Elution: Eluant 2% Glacial Acetic Acid / Methanol
Volume 1.5 ml
Incubation 2 min
Rate 1 ml / min

Analysis: by FPIA



Dr. Mallwitz
Head of Quality Control