


T2-Toxin Standard Solution

known concentration of mycotoxin

Order-No: CH-07-L5-100

Lot: xxx xxx xxx xxx

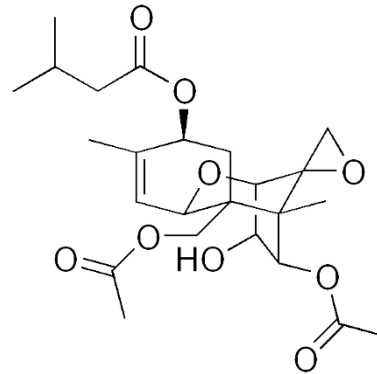

Analyte: T2-Toxin (T2)

Specification:

Substance: T2-Toxin
 Empirical Formula: C₂₄H₃₄O₉
 Concentration: 100 µg/mL
 Diluted in: Acetonitrile
 Molecular Weight: 466,52

CAS-No.: T2-Toxin: 21259-20-1, Acetonitrile: 75-05-8

Volume: 5 mL


Expiry date: 1 year after delivery

Storage conditions: dark, at 2-8°C

Certification: The calibrant is certified on the basis of gravimetric preparation.

Values are based on weight amount, purity and dilution steps, and confirmed by Kinetic Fluorescence Polarization and UV spectroscopy.

Uncertainty ± 1,2 µg/mL in accordance with ISO Guide 31, ISO Guide 35 and Eurachem/CITAG Guides.

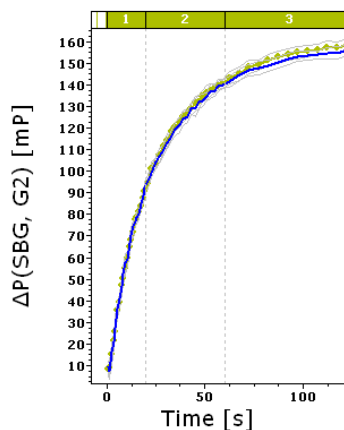

blue: **aokin** T2-Toxin-Standard
 green: external reference standard

Calculation of uncertainty:

(After the concentration of the gravimetric prepared solution was confirmed by kinetic fluorescent polarization, the uncertainty of the calibrant solution was calculated on the basis of preparation)

Uncertainty components	Description	Standard uncertainty (u)	
Purity (P) of solid T2-Toxin	P = 99 % ± 0.5%	u(P) = 0.3%	a
Weighing procedure weighted sample: m _{ws} = 5.1 mg	repeatability: 0.03 mg linearity: 0.012 mg	u(m) = 0.032 mg	b
Dilution procedure volumetric flask: V _f = 50 mL	calibration: 100 mL ± 0.25 mL repeatability: 0.1 mL volume expansion solvent	u(cal) = 0.1 mL u(rep) = 0.1 mL u(Vol.exp.) = 0.119 mL u(V) = 0.12 mL	c d e f

^a Maximum tolerance of purity (rectangular distribution) was divided by $\sqrt{3}$

^b Estimation of this u-value is based upon the values for repeatability and linearity described in the user manual of the microbalance

^c A triangular distribution (division by $\sqrt{6}$) was chosen for the calculation of u(cal)

^d Based on a series of ten fill and weigh experiments on a typical 500 mL flask; the value was used directly as a standard deviation

^e Based on the density of 0.7857 g/cm³ at temperature T = 20°C and a maximum temperature variation of ± 3°C, of volume expansion, relative volume expansion coefficient of acetonitrile is $1370 \cdot 10^{-6}/^{\circ}\text{C}$, volume expansion term (rectangular distribution) was divided by $\sqrt{3}$

^f The three contributions are combined to give the $u(V) = \sqrt{u(\text{cal})^2 + u(\text{rep})^2 + u(\text{Vol.exp.})^2}$

Calculation of the combined uncertainty u_c and the expanded standard uncertainty U:

$$c_{\text{toxin}} = \frac{m_{\text{ws}} \times P}{V_f} = \frac{5.1 \text{ mg} \times 0.99}{50 \text{ ml}} = 0.1 \text{ mg/ml} = 100 \text{ } \mu\text{g/ml}$$

$$\frac{u_c(c_{\text{toxin}})}{c_{\text{toxin}}} = \sqrt{\left[\frac{u(P)}{P}\right]^2 + \left[\frac{u(m)}{m_{\text{ws}}}\right]^2 + \left[\frac{u(V)}{V_f}\right]^2} = \sqrt{\left[\frac{0.3}{99}\right]^2 + \left[\frac{0.032}{5.1}\right]^2 + \left[\frac{0.12}{50}\right]^2} = 0.006$$

$$u_c(c_{\text{toxin}}) = c_{\text{toxin}} \times 0.006 = 100.0 \text{ } \mu\text{g/ml} \times 0.006 = 0.6 \text{ } \mu\text{g/ml}$$

Calculation of expanded standard uncertainty U using a coverage factor k = 2

$$U(c_{\text{toxin}}) = u_c(c_{\text{toxin}}) \times 2 = 0.6 \text{ } \mu\text{g/ml} \times 2 = 1.2 \text{ } \mu\text{g/ml}$$

Discussion of traceability:

This calibrant is certified on the basis of gravimetric preparation. Thus the certified values (mass concentrations of T2-Toxin) is based on the weighed amount of the starting material and are therefore traceable to the stated purity of the solid mycotoxin. High purity material represents a practical realization of concentration units, through conversion of mass to molar quantity.



Danger

H225-H302-H312-H319-H332

P210-P280-P305 + P351 + P338

Danger

Contains: T2-Toxin, Acetonitrile

Volume: 5 mL

Highly flammable liquid and vapour. Harmful if swallowed. Harmful in contact with skin. Causes serious eye irritation. Harmful if inhaled.

Keep away from heat/sparks/open flames/hot surfaces. - No smoking. Wear protective gloves/protective clothing/eye protection/face protection. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

Normal laboratory safety should be observed.

Aokin AG - 13125 Berlin
Tel: +49 (0) 3094892160