





ImmunoAffinity Column for T2/HT2 toxin
for sample preparation

Order no.: IC-C-78-25
Expiry date: XX/XXXX
Lot: 078XXXXXXXXXX071

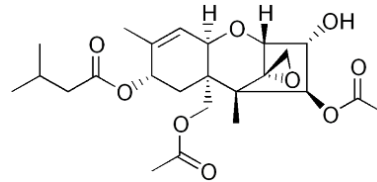

T2-Toxin (T2)

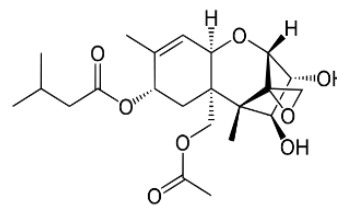

HT2-Toxin (HT2)

Specification:

Recommended loading: < 500 ng T2/HT2
Loading rate: ≤ 1 ml/min
Elution rate: < 1 ml/min

	Recovery of 300 ng
Test results:	T2 ~100 % (n=2)
Specification:	<input checked="" type="checkbox"/> Recovery > 80%

Methodology Testing:

Sample loading: T2 standard 300 ng T2 respectively HT2 standard
Fluid 10 % Methanol / 90 % 10 mM Phosphatbuffer (1:10)
Volume 10 ml
Rate ~ 1 ml / min

Column washing: Fluid 10 mM Phosphatbuffer
Volume 5 ml
Rate > 1 ml / min

Elution: Eluant 100 % Methanol
Volume 2 ml
Rate ~ 1 ml / min
Incubation 3 min

Analysis: by FPIA

Manual

[aokinImmunoClean columns for the quantification of T2 toxin and HT2 toxin](#)